



PRODUCT CATALOG LYMPHOLOGY

2023 / 2024

BAUERFEIND. BY PEOPLE FOR PEOPLE. FOR A BETTER QUALITY OF LIFE.

Bauerfeind AG is one of the world's leading manufacturers of medical aids. With our products, we enable people to be more active and thus make an important contribution toward maintaining and restoring health.

MEDICAL KNOW-HOW AND TECHNICAL EXPERTISE

We have been developing compression products for phlebological and lymphatic conditions for decades, approaching our production and service with untiring dedication and our development with ingenuity.

EFFECTIVE VARIETY

Thanks to our custom-made **VenoTrain curaflow** flat knit products, we can offer patients suffering from lymphedema or lipedema perfectly tailored solutions, with eleven different colors (caramel, navy, anthracite, black, cream, espresso, royal, wild berry, deep petrol, rich olive, and flashy papaya) and more than 1,000 different ways of combining products and additional features. The beneficial therapeutic features of the range include the knitted fabric's high microfiber content, the flat sensitive seams and special functional zones, such as the extra-soft area behind the knee.

A wide range of benefits is combined in our **VenoTrain impuls** circular knit compression stocking: Its breathable microfiber and firm compression thread offer unique wearing comfort with a high degree of effectiveness. Only very little movement is needed for the stocking to exert the required therapeutic pressure.

The circular knit **VenoTrain soft arm sleeve** is used to treat mild arm lymphedema and to prevent imminent arm lymphedema after removal of the lymph node in the armpit or in cases of radiation treatment.*

* For a complete list of all indications, please see the included VenoTrain instructions for use.

UNCOMPROMISING QUALITY

Bauerfeind products don't simply meet standards – they fulfill the most stringent requirements when it comes to quality, effectiveness and wearing comfort.

This is confirmed by inspection procedures for which our products are permitted to bear the appropriate seals and marks after successful testing.

All medical VenoTrain compression stockings are CE-certified, thus complying with current applicable statutory regulations. They also bear the RAL quality label as well as the OEKO-TEX STANDARD 100 quality seal.


From the thread to the completed product, our medical compression stockings are developed and manufactured in Germany.





THE BAUERFEIND SERVICE

RELIABLE SUPPORT AND QUALITY

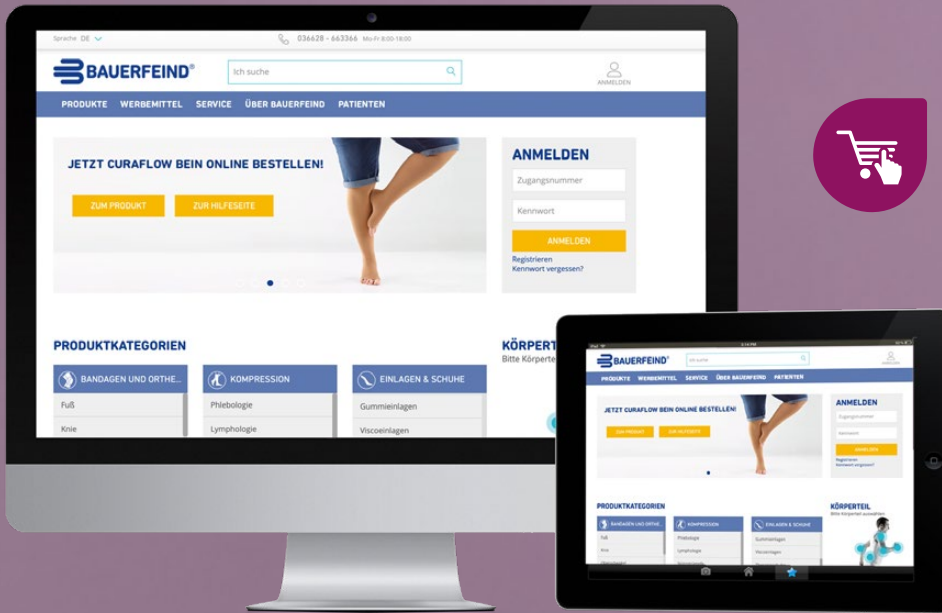

Order online now:
experts.bauerfeind.com

SERVICE GUARANTEES FOR MEDICAL RETAILERS IN GERMANY



FREE TELEPHONE AND FAX HOTLINES

We are happy to help you!



100% SERVICE AVAILABILITY*

As competent contact partners, we are available by phone from 8.00 a.m. to 6.00 p.m. from Monday to Friday. You can reach us around the clock by fax, by e-mail, or via electronic order systems.

* Except in the event of system faults and maintenance work.



ELECTRONIC ORDER SYSTEMS

Place your orders easily, quickly and with reduced shipping costs: Using wheel-it or the various functions available in our online shop, such as the order status and availability display.



SHIPPING ON THE SAME BUSINESS DAY

If you place your order by 3.00 p.m., we will ship in-stock items on the same day (large orders excluded). We are happy to accommodate your express delivery requirements, too.



EXPRESS MANUFACTURE FOR THE QUICK SUPPLY OF CUSTOM-MADE PRODUCTS

Urgent orders for custom-made VenoTrain impuls products (caramel) that are placed before 10:00 a.m. will be completed and shipped on the same business day.



NO SHIPPING COSTS

We do not charge any shipping costs for order values over €200 or for custom-made solutions. You can include all items from any of the product ranges in a single order.



PRECISE DELIVERY DATE

We will specify a precise delivery date for your orders.



FREE GOODS COLLECTION

If you fax us a completed return form, we will come and collect the goods to be returned (subject to inspection) – at our own expense.



365 DAY RETURNS

In accordance with our General Terms and Conditions, we will accept returns of stock items in a faultless condition within 365 days following delivery.



RETURNS PROCESSED WITHIN FIVE DAYS

We will process your returns within five working days of receiving the goods.

CONTENTS

BAUERFEIND ACADEMY
COMPETENT. GOAL-ORIENTED.
SUCCESSFUL.

Page 6

**COLOR OVERVIEW FOR
LYMPHOLOGY PRODUCTS**
THE RIGHT STOCKING FOR
EVERYONE

Page 8

PRODUCT FEATURES
DESIGNED FOR EVERYDAY LIFE

Page 8

EDGING AND TOP BANDS
PERFECT SUPPORT WITH A
HIGH DEGREE OF COMFORT

Page 9

LEG & FOOT



- 11 Care system
- 12 VenoTrain® impuls
Standard and custom-made products
- 14 VenoTrain® delight
Standard and custom-made products
+ Additional features
- 18 VenoTrain® curaflow
Made-to-measure product
+ Additional features
- 25 Types of gusset

ARM & HAND



- 27 Care system
- 28 VenoTrain® soft
Standard products
- 30 VenoTrain® curaflow
Made-to-measure product
+ Additional features for arm sleeves

ACCESSORIES

- 35 VenoTrain® Skin care
- 36 Accessories
- 37 VenoTrain® glider

PRODUCT SELECTION AND ORDER

- 38 Alterations | Repairs
- 40 Order pads
- 41 Additional Bauerfeind treatment areas

BAUERFEIND ACADEMY

COMPETENT. GOAL-ORIENTED. SUCCESSFUL.

Would you like customers to clearly favor your medical retail store and reward your team's expertise?

We will be happy to support you with individual training sessions that include many practical examples.



OUR TRAINING OFFER FOR YOU

Basic and advanced trainings in all business units for example:

- Basic lymphology training
- Basis phlebology training

Take a look at
our website and
see our new
E-Learnings.





BENEFITS AT A GLANCE

+ TRAINING WITH STRONG PRACTICAL RELEVANCE

The training sessions offered by our academy will always include medical information to improve understanding and ensure increased technical expertise.

Specifically and practically educated trainers, such as orthotists, physical therapists, and lymph therapists, will provide professional guidance through the different subjects, while always considering the relevant target group.

+ ONLINE TRAINING

Would you like to learn in a time-saving and cost-efficient way?










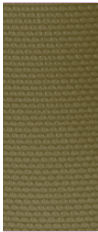

Choose a suitable training session from our wide range of online options.

If you cannot see the subject you are interested in, we will be happy to create a tailor-made online training session for you and to provide an individual quotation.

For further information and registration forms, please contact your local Sales Representative.

COLOR OVERVIEW FOR LYMPHOLOGY PRODUCTS

THE RIGHT PRODUCT FOR EVERYONE






VenoTrain®	STANDARD COLORS				ADDITIONAL COLORS			TREND COLORS NEW			
	Caramel	Navy	Anthracite	Black	Cream	Espresso	Royal	Wild Berry	Deep Petrol	Rich Olive	Flashy Papaya
											
impuls	●	●	●	●				●	●		
delight	●			●							
curaflow	●	●	●	●	●	●	●	●	●	●	●
soft arm	●			●							



PRODUCT FEATURES

DESIGNED FOR EVERYDAY LIFE

VenoTrain curaflow

-  Smooth knitted fabric makes the stockings easy to put on and comfortable to wear
-  Extra-soft popliteal area and inner elbow for uninhibited mobility
-  Smart longitudinal structure for perfect lymph transport
-  Constant pressure and an ideal pressure profile with a tangible effect throughout the time worn, pleasantly soft, fine double-stretch knit
-  Perfect fit for individual body dimensions thanks to anatomically tailored stitched products

EDGING AND TOP BANDS

PERFECT SUPPORT WITH A HIGH LEVEL OF COMFORT

EDGING



COMPRESSIVE EDGING

Edging with compression knit without using silicone or additional seams.

Available for
VenoTrain curaflow



LOOPED ELASTIC DRAWBAND

This edging is only available for pantyhose and is a standard feature.

Available for
VenoTrain curaflow
and VenoTrain impuls



WAIST BAND

Flat edging with an elastic band without silicone coating that offers exceptional wearing comfort.

Available for
VenoTrain curaflow



POROUS EDGING

Reduced pressure at the beginning and/or the end of the stocking; smooth transition when several products overlap.

Available for
VenoTrain delight and
VenoTrain curaflow

SILICONE TOP BANDS



SILICONE MICRO DOTS

The highly elastic silicone top band with its small, flat silicone dots is very inconspicuous, even when worn under tight clothing. The skin-friendly material with its soft surface adapts perfectly to the leg shape.

Available for VenoTrain curaflow



SILICONE DOTTED BAND

Our classic silicone top band provides patients with secure, skin-friendly support.

Available for VenoTrain impuls,
VenoTrain delight and VenoTrain curaflow



STRIPED SILICONE TOP BAND

This silicone top band is coated with three flat stripes and provides stronger support for firmer tissue and larger circumferences.

Available for VenoTrain curaflow



CARE SYSTEM^{*}

LEG AND FOOT

The Bauerfeind care system will help you find the right VenoTrain compression stocking for your customers. The selection of the compression stocking and the compression class should always take into account the individual patient situation as well as the patient's specific needs. If in doubt, please contact us.

LYMPHEDEMA

Stages

Findings / symptoms

PRODUCT RECOMMENDATION

VenoTrain[®] impuls
circular knit

VenoTrain[®] delight
flat knit

VenoTrain[®] curaflow
flat knit

STAGE 0 LATENCY STAGE



The lymphatic system is damaged, but there are no visible or palpable edemas yet.

STAGE 1 SPONTA- NEOUSLY REVERSIBLE



The tissue is soft and can be pressed down gently with a visible indentation. The lymphedema subsides when the leg or arm is elevated.

STAGE 2 NOT SPONTA- NEOUSLY REVERSIBLE



The tissue is hardened and it is either difficult or impossible to press down on it and create an indentation. There may be restricted movement.

STAGE 3 NOT SPONTA- NEOUSLY REVERSIBLE



A significant increase in edema size, sometimes with extreme dimensions.

LIPPEDEMA

Stages

Findings / symptoms

PRODUCT RECOMMENDATION

VenoTrain[®] impuls
circular knit

VenoTrain[®] delight
flat knit

VenoTrain[®] curaflow
flat knit

STAGE 1



The surface of the skin is smooth, the subcutaneous adipose tissue is thickened, and the fat structure is finely knotted.

STAGE 2



The surface of the skin is uneven and the fatty structure features coarse nodules.

STAGE 3



In addition to Stage 2 symptoms, the tissue is coarser and harder as well. There are also lobes of fat causing extensive deformation.

Ccl 1: 18–21 mmHg

- Light compression
- Treatment of childhood lymphedema, edema caused by inactivity

Ccl 2: 23–32 mmHg

- Moderate compression
- Moderate surface effect
- Treatment of primary and secondary lymphedema of the leg (stages 1 and 2), lipedema, lipedema with secondary lymphedema

Ccl 3: 34–46 mmHg

- Firm compression level
- Significant surface effect
- Treatment of primary and secondary lymphedema of the leg (stages 2 and 3), lipedema, lipedema with secondary lymphedema

Ccl 4: from 49 mmHG

- Extra-firm compression level
- Treatment of primary and secondary lymphedema of the leg (stages 2 and 3), lipedema, lipedema with secondary lymphedema

* Graphic representation

VenoTrain® impuls

Highly effective for early-stage lymphedema



The model is wearing
thigh-high stockings (AG)
in navy

Care instructions



V6+ STANDARD SIZE CHART

	CIRCUM- FERENCE	I		II		III		IV		V		VI	
		REGULAR	PLUS	REGULAR	PLUS	REGULAR	PLUS	REGULAR	PLUS	REGULAR	PLUS	REGULAR	PLUS
Length IG	cG	45–53	51–59	48–56	54–62	51–59	57–66	54–63	60–69	57–66	63–73	60–70	66–76
short 62–70	cF	39–45	43–50	42–49	46–53	45–52	49–57	48–56	52–60	51–59	55–64	54–63	58–67
long 70–80	cE	32–37	35–40	34–39	37–42	36–41	39–44	38–44	41–47	40–46	43–49	42–48	45–51
Type AG	cD	28–32	30–34	30–34	32–37	32–37	34–39	34–39	36–41	36–41	38–43	38–43	40–46
Length ID	cC	29–33	31–35	31–35	33–38	33–38	35–40	35–40	37–42	37–42	39–44	39–44	41–47
short 35–40	cB1	23–27	24–28	25–29	26–30	27–31	28–32	29–34	30–35	31–36	32–37	33–38	34–39
long 40–44	cB	19–20.5	20.5–22	20.5–22	22–24	22–24	24–26	24–26	26–28.5	26–28.5	28.5–31	28.5–31	
Type AD	cY	25–29	26–31	26–31	28–33	28–33	30–35	30–35	32–38	32–38	35–41	35–41	
	cA	19–22	20–23	20–23	22–25	22–25	24–27	24–27	26–30	26–30	28–33	28–33	


* All measurements in cm



Lower leg stocking (AD) in navy with closed toe



Thigh-high stocking (AG) in caramel with open toe

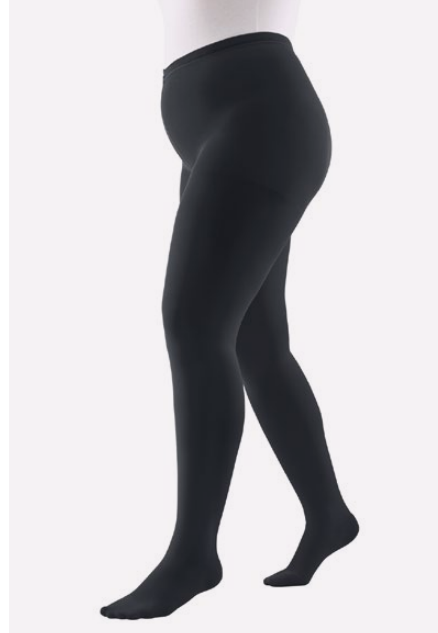

Pantyhose (AT) in anthracite with closed toe

FEATURES

- The short-stretch elastic knitted fabric generates a high working pressure with even the slightest muscle movements.
- Skin-friendly microfiber material makes the product breathable and ensures climate control
- Highly differentiated V6+ sizing system (I – VI)

ORDERING OPTIONS

Sizing system: V6+ (I – VI)

Compression classes:

Ccl 2 (23–32 mmHg)

Ccl 3 (34–46 mmHg)

Toe: open, closed

Additional features: see price list

MATERIAL CONTENT

Ccl 2: 80% polyamide, 20% elastane

Ccl 3: 75% polyamide, 25% elastane

Colors:

Standard colors

CARAMEL



NAVY



ANTHRACITE



BLACK



Trend colors

WILD BERRY



DEEP PETROL



● STANDARD SIZES

○ CUSTOM-MADE SOLUTIONS

TOP BAND	AD	AG	AGHB	AT	ATH	ATX
EDGING						
Top band	● ○			○	○	○
Silicone dots	● ○	● ○				
Garter belt			● ○			

VenoTrain® delight

For intensive therapy



The model is wearing
thigh-high stockings (AG)
in caramel

Care instructions



STANDARD SIZE CHART

		CIRCUMFERENCE	I	II	III	IV	V	VI
Length /D short 35–39 long 40–44 Type AD		cD	28–34	30–37	32–39	34–41	36–43	38–46
		cC	29–35	31–38	33–40	35–42	37–44	39–47
		cB1	23–28	25–30	27–32	29–35	31–37	33–39
		cB	19–20.5	20.5–22	22–24	24–26	26–28.5	28.5–31
		cY	25–29	26–31	28–33	30–35	32–38	35–41
		cA	19–22	20–23	22–25	24–27	26–30	28–33

* All measurements in cm



Foot straight closed



Silicone dotted band (AD) in caramel

FEATURES

- Deep action for the venous system
- Breathable double-stretch fabric
- Special gripRelief® structure at the top ensures the stocking remains comfortably and securely in place
- Extra-flat seams

ORDERING OPTIONS

Sizing system: Standard and custom-made solutions

Compression classes:

Ccl 2 (23–32 mmHg)

Ccl 3 (34–46 mmHg)

Toe: open, closed, diagonal cut

Additional features: see the following pages and price list

Colors:

Standard colors

CARAMEL



BLACK













MATERIAL CONTENT

Ccl 2: 80% polyamide, 20% elastane

Ccl 3: 79% polyamide, 21% elastane

● STANDARD SIZES

○ CUSTOM-MADE SOLUTIONS

										
TOP BAND EDGING	AD	AF	AG	AT	ET	CT	BD	BG	CG	DG
Top band	●○	○	○	○	○	○	○	○	○	○
Silicone dots (3 cm)	●○									
Silicone dots (5 cm)	●○	○	○					○	○	○
Garter belt			○							
TOE										
open	●○	○	○	○						
closed	●○	○	○	○						
diagonal cut	●○	○	○	○						
ADDITIONAL FEATURE										
T-heel	●○*	○*	○*	○						
Viscoelastic pads	○*	○*	○*							

* Optionally with T-heel or pads

experts.bauerfeind.com/venotrain-delight

ADDITIONAL FEATURES FOR LEG AND FOOT

(VenoTrain delight)

ADDITIONAL FEATURES	APPLICATIONS	MODES OF ACTION
Garter belt 	<ul style="list-style-type: none"> • Intolerability to silicone top band 	<ul style="list-style-type: none"> • Support function • Ensures that the product always remains in the correct position
Slanted foot 	<ul style="list-style-type: none"> • Edematous dorsum of the foot • Toe misalignments • Skin and nail conditions in the toe area 	<ul style="list-style-type: none"> • Required pressure gradient at the metatarsophalangeal joints • Anatomical fit • Distal fixation • Prevention of circular constriction and development of forefoot edema
Straight foot (closed) 	<ul style="list-style-type: none"> • Edematous dorsum of the foot • Toe and foot misalignments • Narrow foot tapering toward the heel 	<ul style="list-style-type: none"> • Stops the foot knit slipping backward • Prevents backlog in the toes • Very well tolerated soft toe
T-heel 	<ul style="list-style-type: none"> • Skin folds at Y • Pronounced anatomical heel • Predominantly sedentary tasks 	<ul style="list-style-type: none"> • No creasing at the heel or instep • Improved therapeutic gradient in the ankle area • Better joint functionality
Pad 	<ul style="list-style-type: none"> • Edema and hardening in the ankle area • Depending on indication 	<ul style="list-style-type: none"> • Targeted pressure increase • Prevention of edemas and hardening • Loosening fibroses



VenoTrain[®] curaflow

Compression stocking for treating
lymphedema and lipedema of the leg

Training video

video.bauerfeind.com



The model is wearing
thigh-high stockings (AG) in
flashy papaya with a silicone
dotted band

Care instructions





Thigh-high stocking (AG) in caramel with open toe



Pantyhose (AT) in caramel with waist band

FEATURES

- Smooth microfibers make the product convenient to put on
- Skin-friendly microfiber material makes the product breathable and ensures climate control
- Flat sensitive stitching
- Extra-soft toe
- Anatomical, extra-soft area behind the knee for unrestricted freedom of movement
- Custom-made solutions to suit individual requirements

ORDERING OPTIONS

Sizing system: Custom-made

Compression classes:

Ccl 1 (18–21 mmHg) to
Ccl 4 (from 49 mmHg)

Toe: open, diagonal cut, closed, closed at an angle

Additional features: see the following pages and price list

MATERIAL CONTENT

Ccl 1–2: 73% polyamide, 27% elastane

Ccl 3: 70% polyamide, 30% elastane

Ccl 4: 65% polyamide, 35% elastane

Seamless toe cap:

Ccl 1: 87% polyamide, 13% elastane

Ccl 2: 85% polyamide, 15% elastane

Ccl 3: 86% polyamide, 14% elastane

Toe cap with seam:

Ccl 1–2: 73% polyamide, 27% elastane

Ccl 3: 70% polyamide, 30% elastane

Colors:

Standard colors

CARAMEL



NAVY



ANTHRACITE



BLACK



Additional colors

CREAM



ESPRESSO



ROYAL



Trend colors

WILD BERRY



DEEP PETROL




RICH OLIVE



FLASHY PAPAYA



CUSTOM-MADE SOLUTIONS

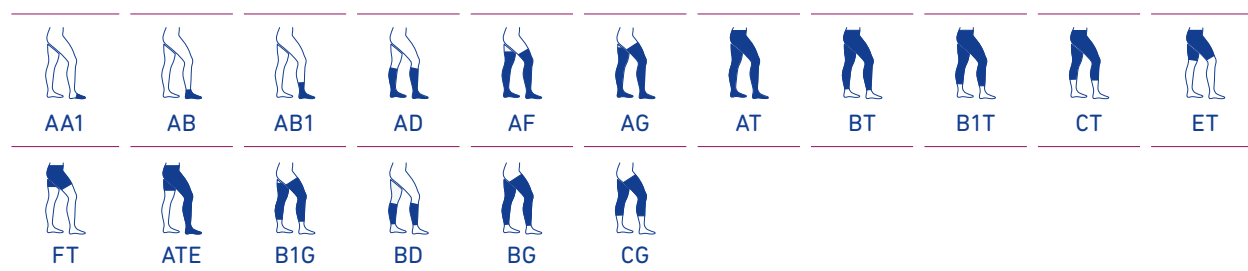

<https://experts.bauerfeind.com/curaflow-leg>

Stage 0 lymphedema

Lymphedema / lipedema stage 1

Lymphedema / lipedema stage 2

Lymphedema / lipedema stage 3

ADDITIONAL FEATURES FOR LEG AND FOOT

(VenoTrain curaflow)

Thanks to the custom-made solutions available, VenoTrain curaflow products can be individually adapted to suit each wearer's requirements. The product features shown here are the most popular additional options.




Our extensive range of custom-made solutions offers a total of over 1,000 product combination options. Please refer to our current price list for more additional features.

ADDITIONAL FEATURES	APPLICATIONS	MODES OF ACTION
Straight foot (open) 	<ul style="list-style-type: none"> • In cases of inconspicuous anatomy 	<ul style="list-style-type: none"> • For improved air circulation in the forefoot area
Slanted foot (open) 	<ul style="list-style-type: none"> • Edematous dorsum of the foot • Toe misalignments • Skin and nail conditions in the toe area 	<ul style="list-style-type: none"> • Required pressure gradient at the metatarsophalangeal joints • Anatomical fit • Distal fixation • Prevention of circular constriction and development of forefoot edema
Straight foot (closed) 	<ul style="list-style-type: none"> • Edematous dorsum of the foot with a relatively straight toe line 	<ul style="list-style-type: none"> • Stops the foot knit slipping backward • Prevents backlog in the toes • Very well tolerated soft toe
Slanted foot (closed) 	<ul style="list-style-type: none"> • Edematous dorsum of the foot • Toe and foot misalignments 	<ul style="list-style-type: none"> • Required pressure gradient at the metatarsophalangeal joints • Anatomical fit • Stops the foot knit slipping backward • Prevents backlog in the toes • Very well tolerated soft toe

ADDITIONAL FEATURES	APPLICATIONS	MODES OF ACTION
Hallux 	<ul style="list-style-type: none"> • Bunion (hallux valgus) 	<ul style="list-style-type: none"> • Perfect fit in the area of the big toe • Pressure relief
90° heel 	<ul style="list-style-type: none"> • Skin folds at Y • Pronounced anatomical heel • For patients in wheelchairs and those who spend much of their time sitting down • For performing everyday activities 	<ul style="list-style-type: none"> • No creasing at the heel or instep • Improved therapeutic gradient in the ankle area • Better joint functionality
Instep widening 	<ul style="list-style-type: none"> • Overhang at Y • High instep • Lymphedema stage end of Stage 2 / beginning of Stage 3 	<ul style="list-style-type: none"> • Reduced pressure on the dorsiflexor tendon • Improved compression gradient along the dorsum of the foot
Seamless toe cap 	<ul style="list-style-type: none"> • Forefoot and toe lymphedema • Edema of the dorsum of the foot 	<ul style="list-style-type: none"> • Prevention of edema in the forefoot and toe area • Improved compliance, thanks to wearing comfort
Toe cap 	<ul style="list-style-type: none"> • Forefoot and toe lymphedema • Edema of the dorsum of the foot 	<ul style="list-style-type: none"> • Prevention of edema in the forefoot and toe area
Toe cap with rings 	<ul style="list-style-type: none"> • Forefoot and toe lymphedema • Edema of the dorsum of the foot 	<ul style="list-style-type: none"> • Prevention of edema in the forefoot and toe area

ADDITIONAL FEATURES	APPLICATIONS	MODES OF ACTION
Toes sewn onto the stocking 	<ul style="list-style-type: none"> • Forefoot and toe lymphedema • Edema of the dorsum of the foot 	<ul style="list-style-type: none"> • Prevention of edema in the forefoot and toe area • No overlapping of the stocking and toe section
Integrated pad 	<ul style="list-style-type: none"> • Edema and hardening in the ankle area 	<ul style="list-style-type: none"> • Targeted pressure increase • Prevention of edemas and hardening
Integrated lymph pad 	<ul style="list-style-type: none"> • Locally aggravated edema or hardening • Anatomical deviations from standard shapes 	<ul style="list-style-type: none"> • Targeted pressure increase • Prevention of edemas and hardening
Lining 	<ul style="list-style-type: none"> • Skin irritation in bending areas, such as the ankle and knee, parchment skin 	<ul style="list-style-type: none"> • Protection of sensitive skin from excessive strain • Reduction of rubbing and mechanical load • Securing of pads
Functional zone Standard feature on the knee for all stockings 	<ul style="list-style-type: none"> • Generally for lymphedema and lipedema patients 	<ul style="list-style-type: none"> • Extra-soft area behind the knee to reduce creasing • Anatomical fit for improved joint functionality
Slanted edging on the knee 	<ul style="list-style-type: none"> • Generally for lymphedema and lipedema patients • Tapered anatomical shape • Soft skin tissue 	<ul style="list-style-type: none"> • Improved therapeutic gradient at the ends of the product (AD) • Anatomical fit

ADDITIONAL FEATURES	APPLICATIONS	MODES OF ACTION
Slanted edging on the thigh	<ul style="list-style-type: none">• Generally for lymphedema and lipedema patients• Tapered anatomical shape• Soft skin tissue	<ul style="list-style-type: none">• Prevents circular constriction of lymphatic vessels• Improved therapeutic gradient at the ends of the product (AG)• Anatomical fit
Leg elevation at the front only available in combination with slanted edging	<ul style="list-style-type: none">• Patients with voluminous thigh marked by edema	<ul style="list-style-type: none">• Prevents constriction• Physiological pressure gradient up to the groin
Silicone top band section	<ul style="list-style-type: none">• Pronounced tapered leg shape• Very soft tissue	<ul style="list-style-type: none">• Extra grip
Men's gusset Standard	<ul style="list-style-type: none">• Pantyhose for men	<ul style="list-style-type: none">• Improved fit in the panty area

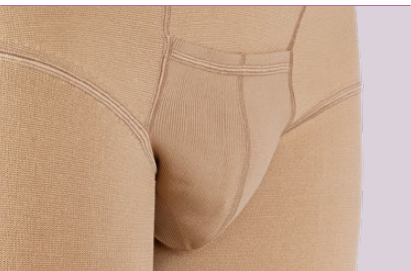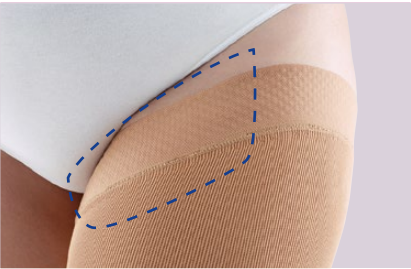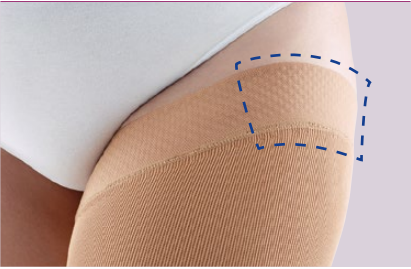

ADDITIONAL FEATURES	APPLICATIONS	MODES OF ACTION
Stoma opening 	<ul style="list-style-type: none"> • Abdomen 	<ul style="list-style-type: none"> • Allows access to external drainage / supply
ATU gusset 	<ul style="list-style-type: none"> • Pregnancy 	<ul style="list-style-type: none"> • Improved fit in the panty area • Reduced compression in the abdominal region
Garter belt 	<ul style="list-style-type: none"> • Intolerability to silicone top band 	<ul style="list-style-type: none"> • Holds the compression stocking in place
Zip fastening 	<ul style="list-style-type: none"> • For major differences between hips and waist circumferences • In cases of weak paresis 	<ul style="list-style-type: none"> • Makes donning and doffing of the product easier

TYPES OF GUSSET

Standard gusset without transverse opening




This is the standard product which, with its single fabric layer, is suitable for the use of hygiene products and prevents constriction in the genital area.

Standard gusset with transverse opening



Thanks to the double fabric layer, rubbing or irritation in the genital area will be prevented. Furthermore, it helps absorb sweat and discharge, thus reducing the risk of infection.

Compressive gusset without transverse opening

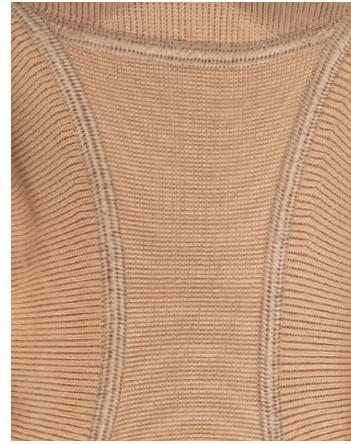

The excellent fit prevents slipping and ensures the knitted fabric remains securely in position.

Compressive gusset with transverse opening



The excellent fit prevents slipping and ensures the knitted fabric remains securely in position. The double fabric layer also prevents rubbing and skin irritation in the genital area.

Overlapping gusset



Allows post-operative use in patients with catheters and makes using the restroom easier.



CARE SYSTEM*

ARM AND HAND

The Bauerfeind care system will help you find the right VenoTrain compression stocking for your customers. The selection of the compression stocking and the compression class should always take into account the individual patient situation as well as the patient's specific needs. If in doubt, please contact us.

LYMPHEDEMA

Stages

Findings / symptoms

PRODUCT RECOMMENDATION

VenoTrain® soft
circular knit

VenoTrain® curaflow
flat knit



STAGE 0 LATENCY STAGE

The lymphatic system is damaged, but there are no visible or palpable edemas yet.



STAGE 1 SPONTANEOUSLY REVERSIBLE

The tissue is soft and can be pressed down gently with a visible indentation.



STAGE 2 NOT SPONTANEOUSLY REVERSIBLE

The tissue is hardened and it is either difficult or impossible to press down on it and create an indentation.



STAGE 3 NOT SPONTANEOUSLY REVERSIBLE

A significant increase in edema size, sometimes with extreme dimensions.

LIPDEMA

Stages

Findings / symptoms

PRODUCT RECOMMENDATION

VenoTrain® soft
circular knit

VenoTrain® curaflow
flat knit



STAGE 1

The surface of the skin is smooth, the subcutaneous adipose tissue is thickened, and the fat structure is finely knotted.



STAGE 2

The surface of the skin is uneven and the fatty structure features coarse nodules.



STAGE 3

In addition to Stage 2 symptoms, the tissue is coarser and harder as well. There are also lobes of fat causing extensive deformation.

COMPRESSION CLASSES

Ccl 1: 18 – 21 mmHg

- Light compression

Ccl 2: 23 – 32 mmHg

- Moderate surface effect
- Treatment of primary and secondary lymphedema of the arm (stages 1 and 2), lipedema, lipedema with secondary lymphedema

Ccl 3: 34 – 46 mmHg

- Significant surface effect
- Treatment of primary and secondary lymphedema of the arm (stages 2 and 3), lipedema, lipedema with secondary lymphedema

* Graphic representation

VenoTrain® soft arm sleeve

Freedom of movement during everyday activities




STANDARD SIZE CHART VenoTrain® soft arm

Length /D

short 39–44

long 44,1–50

Type CG

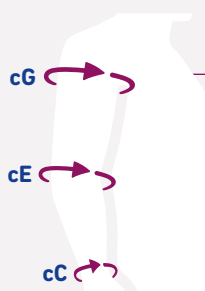

CIRCUM- FERENCE	S MIN. – MAX.	M MIN. – MAX.	L MIN. – MAX.	XL MIN. – MAX.
cG	24 – 30	27 – 33	29 – 36	32 – 40
cE	22 – 26	24 – 29	27 – 32	30 – 35
cC	14 – 15	15 – 16.5	16.5 – 18	18 – 20

* All measurements in cm



Functional zone at the crook of the elbow in beige



Top band in beige

FEATURES

- Comfortable inner elbow zone for greater freedom of movement, wearing comfort and reduced feeling of pressure
- Noticeable effect throughout the entire wearing period
- Gentle on the skin, soft, and breathable

ORDERING OPTIONS

Sizing system:

VenoTrain soft arm (S–XL)

Compression classes:

Ccl 1 (18–21 mmHg)

Ccl 2 (23–32 mmHg)

MATERIAL CONTENT

Ccl 1: 74% polyamide, 26% elastane

Ccl 2: 69% polyamide, 31% elastane

Colors:

Standard colors

BEIGE



CARAMEL



BLACK



STANDARD SIZES

TOP BAND EDGING



CG

Top band



Silicone dots



VenoTrain[®] curaflow

Compression stocking for treating lymphedema and lipedema of the arm and hand



Training video
video.bauerfeind.com



The model is wearing an arm sleeve (CG) and the seamless glove (AC1) in royal

Care instructions





Seamless glove (AC1) in royal with open fingers



Arm sleeve (CG) in caramel

FEATURES

- Smooth microfibers make the product convenient to put on
- Skin-friendly microfiber material makes the product breathable and ensures climate control
- Anatomically contoured, extra-soft inner elbow area for greater freedom of movement, wearing comfort and reduced feeling of pressure
- Custom-made solutions to suit individual requirements

ORDERING OPTIONS

Sizing system: Custom-made

Compression classes:

Ccl 1 (18–21 mmHg) to

Ccl 3 (34–46 mmHg)

Additional features: see the following pages and price list

MATERIAL CONTENT

Arm:

Ccl 1: 77% polyamide, 23% elastane

Ccl 2: 77% polyamide, 23% elastane

Ccl 3: 73% polyamide, 27% elastane

Seamless hand:

Ccl 1: 87% polyamide, 13% elastane

Ccl 2: 85% polyamide, 15% elastane

Ccl 3: 86% polyamide, 14% elastane

Hand with seam:

Ccl 1: 73% polyamide, 27% elastane

Ccl 2: 73% polyamide, 27% elastane

Ccl 3: 70% polyamide, 30% elastane

Colors:

Standard colors

CARAMEL



NAVY



ANTHRACITE



BLACK



Additional colors

CREAM



ESPRESSO



ROYAL



Trend colors

WILD BERRY



DEEP PETROL




RICH OLIVE



FLASHY PAPAYA






CUSTOM-MADE SOLUTIONS

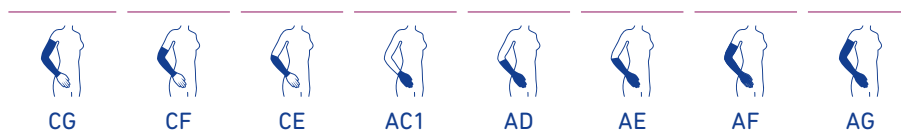

ADDITIONAL FEATURES FOR ARM AND HAND

Thanks to our portfolio of custom-made solutions, VenoTrain curaflow products can be individually adapted to suit each wearer's requirements. The product features shown here are the most popular additional options. Our extensive range of custom-made solutions offers a total of over 1,000 product combination options.

Please refer to our current price list to find out more about additional features or order them from the Bauerfeind online shop.

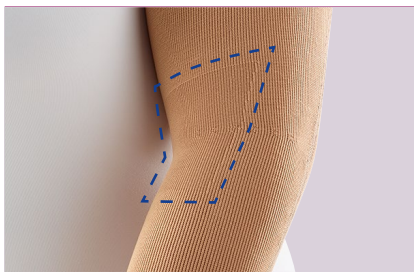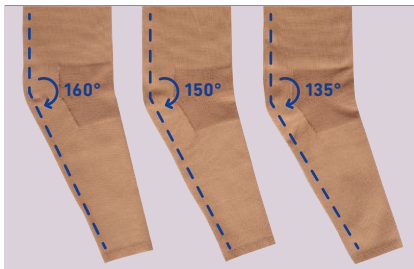ADDITIONAL FEATURES	APPLICATIONS	MODES OF ACTION
<p>Functional zone Standard feature on all arm sleeves</p> 	<ul style="list-style-type: none"> • Generally for lymphedema and lipedema patients 	<ul style="list-style-type: none"> • Integrated functional zone to reduce creasing • Anatomical fit for improved joint functionality
<p>Elbow 160°, 150°, or 135°</p> 	<p>160°: as a standard product 150°: ongoing work with a bent elbow 135°: patients who can no longer fully extend the elbow</p>	<ul style="list-style-type: none"> • Prevents creasing in the bending area, thus preventing constriction
<p>Slanted edging</p> 	<ul style="list-style-type: none"> • Generally for lymphedema and lipedema patients • Anatomical shape • Soft skin tissue 	<ul style="list-style-type: none"> • Improved therapeutic gradient at the ends of the product • Anatomical fit • Prevents circular constriction

ADDITIONAL FEATURES	APPLICATIONS	MODES OF ACTION
Silicone top band section 	<ul style="list-style-type: none"> • For very soft tissue • Pronounced tapered arm shape 	<ul style="list-style-type: none"> • Extra grip • Prevents slipping
Zip fastening 	<ul style="list-style-type: none"> • Secondary conditions such as rheumatism, arthritis, etc. • In cases of weak paresis 	<ul style="list-style-type: none"> • Makes donning and doffing of the product easier
Integrated lymph pad 	<ul style="list-style-type: none"> • Locally aggravated edema or hardening • Anatomical deviations from standard shapes 	<ul style="list-style-type: none"> • Targeted pressure increase • Prevention of edemas and hardening
Seamless glove 	<ul style="list-style-type: none"> • Mild to moderate edema of the fingers • Mild to moderate edema on the back of the hand 	<ul style="list-style-type: none"> • Prevention of edemas in fingers and on the back of the hand
Glove with seam 	<ul style="list-style-type: none"> • Moderate to pronounced edema of the fingers • Moderate to pronounced edema on the back of the hand 	<ul style="list-style-type: none"> • Prevention of edemas in fingers and on the back of the hand

VenoTrain® ACCESSORIES

THE PERFECT ADDITION
FOR EVERYDAY USE



VenoTrain® SKIN CARE

SPA TREATMENT FOR THE SKIN

Looking after their skin should now be a priority for your customers because the skin is generally more sensitive in cases of lymphedema or lipedema. Skin care creams prevent the skin from drying out and provide the necessary protection. To ensure the skin stays smooth and beautiful, we have developed two high-quality care products that are perfectly adapted to skin subjected to strain from edema and compression.



VenoTrain® LOTION*

A beneficial treat in the evenings

After a long day of wearing compression garments, the skin deserves a little spa treatment. The care lotion we have developed together with Sebapharma is the ideal choice for this. Its special chamomile extract soothes the skin and reduces skin irritation, thus ensuring perfectly relaxed legs and arms.



BAUERFEIND CARE FOAM*

Quick care in the morning

Gentle care foam that protects and soothes the skin using panthenol. The cream absorbs quickly and can be used before putting on compression garments.

* Not available in all countries

VenoTrain® ACCESSORIES

The handy extras from the VenoTrain accessories range are useful for ensuring that patients feel comfortable using and wearing their VenoTrain compression stockings. They support treatment, make it easier to handle the products and care for the compression stockings as well as your customers' skin.



01



02



03

01 Measurement bag

to help with measurements: tailor's chalk, finger measuring tape, eyeliner pencil, pencil sharpener, tape measure, special gloves for the easy donning of compression stockings, VenoTrain glider, pen

02 Lingerie bag

to protect the compression stockings during every machine wash

03 Bauerfeind special gloves

help with putting on the compression stockings correctly

04 Starter kit

the useful accessories kit provides support during daily wearing of compression stockings (e.g. gloves, lingerie bag)

05 Detergent*

cleans gently, ensuring careful treatment of the products and providing lasting protection



04



05

* Not available in all countries

VenoTrain® glider

COMPACT DONNING AIDS

The VenoTrain glider is a practical donning aid for medical compression stockings.

The lightweight, hard-wearing material makes the stocking glide effortlessly over the arm or leg. It is suitable for compression stockings with an open or closed toe.

Saving space and time, it is the ideal companion when traveling.

FEATURES

- Extremely lightweight, robust, and tear-resistant
- Easy to handle
- Ideal for travel

THE VIDEOS WITH DONNING INSTRUCTIONS CAN BE FOUND AT:

video.bauerfeind.com





ADDITIONAL SERVICES ALTERATIONS | REPAIRS INDIVIDUAL CUSTOMER SERVICE

ALTERATIONS | REPAIRS

● STANDARD SIZES

		soft Arm	impuls Leg
Top of stocking	Replace silicone top band	●	●
	Change top band to silicone top band	●	●
	Silicone top band	●	●
	Remove hip section and sew on silicone top band		●
AGHB	Replace hip strap		●
Shorten stocking	AD / CG by up to 3 cm	●	●
	AG by up to 5 cm		●
Repairs, if possible	Hole	●	●
	Run	●	●
	Seam	●	
	Loose thread	●	●

ORDERING NEW PRODUCTS

○ CUSTOM-MADE SOLUTIONS

	Unit	impuls	delight
Knitted toe area Lengthen or shorten	Item	○	
Toe Seam on the outside	Item	○	○
Toe Make wider, narrower, anatomical	Item	○	○
Pantyhose top Front / back – two seams	Item	○	
AT with men's gusset (without fly or with longitudinal slit)	Item	○	
AT with extra-wide gusset for the fuller figure	Item	○	
Sew in a foam pad	Item		○

ALTERATIONS | REPAIRS

○ CUSTOM-MADE SOLUTIONS

		impuls Leg	delight Leg
Top of stocking	Replace silicone top band	○	○
	Change top band to silicone top band	○	○
	Silicone top band narrower or wider	○	○
	Remove hip section and sew on silicone top band	○	○
AT	Shorten pantyhose top	○	○
	Lengthen pantyhose top (excluding ATX)	○	
	Take in gusset (excluding ATX)	○	
	Widen gusset (excluding ATX)	○	
	Extra-wide gusset for the fuller figure (excluding ATX)	○	
	Alter AT to AG with silicone top band	○	
AGHB	Replace hip strap	○	○
Shorten stocking	AD by up to 3 cm	○	○
	AG by up to 5 cm	○	○
Redye from light to dark*		○	○
Toe alteration	Alter closed toe to open toe	○	○
	Shorten open toe	○	○
Repairs, if possible	Hole	○	○
	Run	○	○
	Seam	○	○
	Loose thread	○	○

ALTERATIONS | REPAIRS

○ VenoTrain® curaflow

Alteration		Treatment	
		Leg	Arm
Pantyhose			
Pantyhose top	Incorporate U part subsequently, for compression pantyhose top only	○	
	Shorten pantyhose top		
	Incorporate stoma opening in pantyhose top		
	Waist edging		
	Gusset		
AT to AG with garter belt		○	
AT to CT, ET or FT		○	
AT to AG with silicone top band		○	
Leg / Arm			
Shorten stocking (chain-stitched edging)		○	○
Silicone top band / waist band / strap	Sew in new sections of silicone top band	○	○
	Sew on new silicone top band / waist band		
	Silicone top band up to 5 cm narrower or wider		
	Alter garter belt to AG with silicone top band		
	Alter shoulder support to CG with silicone top band		
	New strap for garter belt / shoulder support (wear)		
Sew garter belt onto leg stockings up to G		○	
Sew shoulder support onto CG			○
Lymph pad / lining		○	○
Finishing		○	○
Foot (per foot)			
Toe alterations / foot length alteration		○	
Reduce cA dimension		○	
Sew silicone top band pieces in at cA		○	
Detach little toe from toe cap (for toe cap with seam only)		○	
Hand alterations			
Separating fingers from hand (not for seamless AC1)			○
Shorten / take in glove (for gloves with seam only)			○
Incorporate / alter lymph pad subsequently			○
Repairs / adjustments			
Seam repair		○	○
Run repair		○	○
Hole and snag repair		○	○
Redye (for unworn products only)		○	○

* Please do **not** try on stockings – redyeing is no longer possible after contact with the skin.

ORDER PADS

Our special order pads will make your daily work easier.
The different order sheets can also be downloaded at any time at:

bauerfeind.com/downloads

ADDITIONAL BAUERFEIND PRODUCT CATEGORIES



VENOTRAIN MEDICAL COMPRESSION STOCKINGS (circular knit)

are characterized by unique material properties. They improve venous activity and counteract venous leg problems. As medical compression stockings for mild to advanced vein problems, but also as support stockings for work or travel, the VenoTrain product range offers a suitable option. They can easily be combined with various styles, thanks to a wide range of colors and designs.



SUPPORTS AND ORTHOSES

by Bauerfeind help after injuries as well as in cases of acute and chronic orthopedic conditions. They improve the feeling of stability and provide confidence during everyday activities and exercise. The supports and orthoses are made from lightweight and breathable materials to perfectly underpin therapeutic exercise. They are easy to handle during everyday use thanks to innovative solutions such as quick-fit fasteners and integrated donning aids and will win you over with their high-quality design.



A healthy body needs healthy feet. **ORTHOPEDIC FOOT ORTHOSES** simultaneously relieve pressure on the feet, knees, hips and back. They help the foot to function naturally and can therefore positively influence the body's biomechanics. As different feet need different types of support, Bauerfeind offers more than 300 different types of foot orthosis.

For an overview of all treatment areas, please visit:

experts.bauerfeind.com

BAUERFEIND AG TERMS AND CONDITIONS OF SALE AND DELIVERY (GTC)

Prices

The prices quoted are net prices in euros ex works and do not include statutory VAT.

Shipping

Shipping is at the buyer's cost and risk even if the buyer does not bear the transportation cost. Disruptions of operations affecting us or any of our suppliers, lack of suitable personnel, strike, lockout, and force majeure are causes for extending delivery and relieve us of any delivery obligation for the time of the obstruction.

Shipping costs

We calculate our shipping charges as a flat rate of EUR 7.00 for a net total order value of up to EUR 50.00, and EUR 4.00 for a net total order value of up to EUR 200.00. No shipping costs applied, when the order contains custom-made goods. Special shipping requirements will be charged at cost.

Returns

Returns outside the warranty terms are only accepted after prior advice and approval from us, provided the original delivery documents are submitted; it is supposed that the returned goods are in good order. We reserve the right of rejecting the acceptance of any returned consignment. Except where covered by warranty, we charge for processing returns. For returns within one month of the date of delivery, this charge is EUR 10.00 for each return consignment. For goods returned one month or more after the date of delivery, a handling fee of 20 percent of the net order value, or a maximum of EUR 25.00 per product and unit, is charged; for returns six to twelve months after delivery date, the fee is 40 percent of the net order value, or a maximum fee of EUR 50.00 for each product and unit.

Goods returned outside the warranty or 12 months or more after the date of delivery will not be accepted back. Custom made goods are non-returnable.

Payments

If not otherwise agreed individually in advance and in writing, Payments are due immediately following delivery. We offer a 4 percent discount for payments made by direct debit and a 3 percent discount for payments made within 10 days of the invoice date. Otherwise, the period allowed for payment is 30 days net. The buyer shall be considered in default once 30 days have passed since the invoice date. From this point on, any further deliveries will only be made against advance payment. From the time at which payment is delayed, we will charge default interest of eight percent above the applicable base interest rate. Payments by check, bank draft, etc. are deemed to have been made when the amount has finally been credited to our account.

Retention of title

We retain the title to all goods delivered until the buyer has paid all amounts outstanding under our business relations with him. The retention of title shall also apply proportionately to goods which have been processed or otherwise modified. If the goods are sold, all third-party claims are deemed to have been assigned to us without provision up to the total amount owing to us.

Notwithstanding the assignment, the buyer is entitled to the collection of these amounts. Our right to collect the amounts is not affected by the preceding provision; however, we undertake to refrain from collection unless the buyer defaults on payment. If the buyer defaults on payment we can demand that the buyer discloses the assigned claims and the names of the debtors, provides all information required for collection of the debt, delivers all required documents and informs the

debtors (third party) of the assignment. The buyer shall not be entitled to sell, pledge or assign the goods as security. Pledges, confiscation, or any other dispositions taken by a third party shall be promptly communicated to us by the buyer, and all information and documents necessary for us to defend our rights shall be provided. We undertake to release, on the buyer's demand, any security which exceeds a secured claim that has not been satisfied by more than 20 percent.

Complaints

Warranty claims in respect of obvious defects shall be void if we do not receive immediate written notification of these claims in accordance with Section 377 of the German Commercial Code (HGB). Warranty claims shall also lapse one year after the delivery of the purchased goods. If the product is defective, we shall make a replacement delivery or repair it at our discretion. For reasons of hygiene, repairs will only be performed on products that have been washed or are completely hygienic. Upon receipt of unwashed goods, we reserve the right to claim a suitable sum for additional work. The customer may not demand a price reduction or withdraw from the contract provided that we fulfil our duties of subsequent performance, particularly with regard to the repair of defects or the delivery of goods without defect, unless our attempts at subsequent performance have failed. Defects which are caused by normal wear and tear, overuse, misuse, neglect of care, or noncompliance with the instructions for use are not covered by warranty.

Liability

Liability for damages is excluded unless there has been grossly negligent or willful conduct on our part or on the part of our representatives or agents. This limitation of liability shall not apply to damage due to injury to life, body or health, or the violation of a material contractual duty.

Place of jurisdiction

The place of performance for the delivery is our company headquarters. The place of jurisdiction for all disputes arising from the business relation is Gera, Germany.

Severability clause

The event that any individual provision of these terms and conditions of sale and delivery is or becomes invalid, will not affect the validity of the remaining provisions. The invalid provision shall be replaced with a valid provision which reflects the original economic intention of the parties as closely as possible.

Data protection clause / recording of conversations

For quality- and other purposes, we record calls, faxes, e-mails and letters between the buyer and our employees. All recordings serve legal provisions, which have to be met with medical products ordered. Recordings will be deleted automatically in line with legal provisions. Non of the recorded data will be disclosed to third parties, if not required by decision of court or on request of local tax or healthcare authorities. The buyer hereby consents to recordings such.

Please see www.bauerfeind.de/ds for further information on the latest version of Bauerfeind AG's General Data Protection Policy.

Last update: February 2022

A BETTER QUALITY OF LIFE FOR YOUR CUSTOMERS!

BAUERFEIND AG

Triebeser Straße 16
07937 Zeulenroda-Triebes
Germany

P +49 (0) 36628 66-10 00

F +49 (0) 36628 66-19 99

E info@bauerfeind.com

AUSTRIA

Bauerfeind Ges.m.b.H.
Hainburger Straße 33
1030 Wien

P +43 (0) 800 4430-130

F +43 (0) 800 4430-131

E info@bauerfeind.at

BENELUX

Bauerfeind Benelux B.V.
Waarderveldweg 1
2031 BK Haarlem

THE NETHERLANDS

P +31 (0) 23 531-9427

F +31 (0) 23 532-1970

E info@bauerfeind.nl

BELGIUM

P +32 (0) 2 527-4060

F +32 (0) 2 792-5345

E info@bauerfeind.be

BOSNIA AND HERZEGOVINA

Bauerfeind d.o.o.
Meše Selmovića 19
71000 Sarajevo

P +387 (0) 33 710-100

F +387 (0) 33 619-422

E info@bauerfeind.ba

CROATIA

Bauerfeind d.o.o.
Goleška 20
10020 Zagreb

P +385 (0) 1 6542-855

F +385 (0) 1 6542-860

E info@bauerfeind.hr

FRANCE

Bauerfeind France S.A.R.L.
13 rue de la Perdrix
93290 Tremblay-en-France

P +33 (0) 1 4863-2896

F +33 (0) 1 4863-2963

E info@bauerfeind.fr

ITALY

Bauerfeind Italia Srl
Via Cornaggia 58
20092 Cinisello Balsamo (MI)

P +39 02 8977 6310

F +39 02 8977 5900

E info@bauerfeind.it

NORDIC

Bauerfeind Nordic AB
Storgatan 14
114 55 Stockholm

P +46 (0) 774 100 020

E info@bauerfeind.se

REPUBLIC OF NORTH MACEDONIA

Bauerfeind Dooel Skopje
50 Divizija 24 a
1000 Skopje

P +389 (0) 2 3179-002

F +389 (0) 2 3179-004

E info@bauerfeind.mk

SERBIA

Bauerfeind d.o.o.
102 Omladinskih brigada
11070 Novi Beograd

P +381 (0) 11 2287-050

F +381 (0) 11 2287-052

E info@bauerfeind.rs

SINGAPORE

Bauerfeind Singapore Pte Ltd.
Blk 41 Cambridge Road #01-21
Singapore 210041

P +65 6396-3497

F +65 6295-5062

E info@bauerfeind.com.sg

SLOVENIA

Bauerfeind d.o.o.
Dolenjska cesta 242 b
1000 Ljubljana

P +386 (0) 1 4272-941

F +386 (0) 1 4272-951

E info@bauerfeind.si

SPAIN

Bauerfeind Ibérica, S.A.
C / San Vicente Mártir,
nº 71 – 4º – 7ª

46007 Valencia

P +34 96 385-6633

F +34 96 385-6699

E info@bauerfeind.es

SWITZERLAND

Bauerfeind AG
Vorderi Böde 5
5452 Oberrohrdorf

P +41 (0) 56 485-8242

F +41 (0) 56 485-8259

E info@bauerfeind.ch

UNITED ARAB EMIRATES

Bauerfeind Middle East FZ-LLC
Dubai Healthcare City
Building 40, Office 510
Dubai

P +971 4 4335-684

F +971 4 4370-344

E info@bauerfeind.ae

UNITED KINGDOM

Bauerfeind UK
85 Tottenham Court Road
London
W1T 4TQ

P +44 (0) 121 446-5353

F +44 (0) 121 446-5454

E info@bauerfeind.co.uk

USA

Bauerfeind USA, Inc.
75 14th St NE
Suite 2350
Atlanta, GA 30309
P +1 800 423-3405
P +1 404 201-7800
F +1 404 201-7839
E info@bauerfeindusa.com